GBUKBanana
Part of the GBUK Group ●●●●

Repositioning Sling

INSTRUCTIONS FOR USE

	Page
EN Repositioning Sling	4
FR Harnais de repositionnement	6
IT Imbracatura per il riposizionamento	8
ES Cabestrillo de reposicionamiento	10
DE Umlagerungsschlinge	12
DA Repositioneringssslynge	14
NO Reposisjoneringssslynge	16
SE Ompositioneringssele	18
NL Strop voor verplaatsing	20
PL Nosidło do zmiany pozycji	22
PT Cinto de reposicionamento	24
EL Αώρα επανατοποθέτησης	26
TR Yer Değiştirme Askısı	28
CS Zvedací závěs	30
ET Teisaldustropp	32
HU Áthelyező heveder	34
LV Pārvietošanas slings	36
MT Braga għar-Ripożizzjonar	38
RO Ham de repoziționare	40
SL Zanka za premeščanje	42

GBUK Group Ltd.
Woodland House,
Blackwood Hall Business Park,
North Duffield, Selby,
North Yorkshire, YO8 5DD, UK
EU: +44 (0)1757 288 587
US: 18665 289 422
www.gbukbanana.com
info@gbukbanana.com

QualRep Services BV
Utrechtseweg 310 B42
6812 AR Arnhem,
The Netherlands

EN

Repositioning Sling

Instructions for Use

Please read these instructions before use.

INTENDED USE

The Repositioning Sling is intended to be used to hoist a patient in a supine position for transfer, in-bed repositioning or turning.

INDICATIONS FOR USE

The GBUK Banana Repositioning Sling is indicated to be used with patients with limited mobility requiring repositioning, turning or transferring between surfaces.

DESCRIPTION

The GBUK Repositioning sling is made from a polyester fabric with 20 nylon webbing loop attachments located on two inside and two outside channels. The GBUK Repositioning sling has a patient label which contains a loop attachment record which caregivers can use to record different colour loop attachments.

PATIENT TARGET GROUP

This device is suitable to be used for patients who require to be repositioned, turned or transferred between surfaces.

INTENDED USER

GBUK Banana Repositioning Slings are intended to be used by trained healthcare professionals or those with the appropriate training or knowledge.

CONTRAINDICATIONS

Do not use with patients whose weight exceeds the safe working load limit of the repositioning sling or hoist, whichever is lower.

WARNINGS & PRECAUTIONS

- Follow the local manual handling policy from your organization, National Back Exchange, HSE or RCN guidelines.
- Ensure that the correct sling is being used.
- All handlers must be trained and competent when using this device.
- Do not use if packaging is open or damaged.
- Do not use if the wash indicator has turned red, or if the sling shows any sign of damage or wear.
- Do not launder or bleach the device.
- For single patient use only. Cross-patient use poses risk of infection and cross-contamination.
- In the event of a serious incident involving the device the user should report the event to GBUK and the competent authority in which the user/patient is established. In the event of a suspected device fault the user should report the event to GBUK, and the device should be returned to GBUK for examination.
- A risk assessment and visual inspection must be performed by competent person(s) before using the repositioning sling. This should take into account the following risk factors:
 - The physical and medical status of the patient, including skin integrity.
 - The status of the sling, including the appearance of the wash indicator and the presence of any damage, rips, tears, soiling, loose stitching, or any other general signs of wear and tear. Do not use if sling is wet.
 - The compatibility of the repositioning sling with the hoist ensuring that the spreader bar is compatible with the loop attachments
 - The capability of the caregivers using the equipment, and,
 - The situation in which it is being used.
- A risk assessment must be performed by competent person(s) before leaving the repositioning sling in situ.

COMPATIBILITY


GBUK Banana Repositioning sling is suitable for use with a loop spreader bar. Ensure a risk assessment is carried out to ascertain the correct positioning of the loops on the spreader bar.

SAFE WORKING LOAD LIMIT

The GBUK Repositioning Slings have a safe working load of up to 454kg.

WASH INDICATOR

The repositioning sling may be wiped clean if required but should not be laundered. If the sling has been exposed to water, the wash indicator will turn red and the sling must be disposed of.

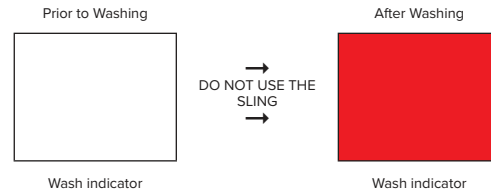OPERATING INSTRUCTIONS:

FITTING THE REPOSITIONING SLING

- Ensure the brakes of the bed/trolley are applied and the bed is of a suitable working height for the caregiver.
- The sling can be inserted by turning the patient on to their side and either folding the sling in half or rolling to the centre point and tucking this under the patients back. Ensuring that the head label is positioned at the head end of the bed, is facing up and is clearly visible.
- The patient can then be turned on to their other side and the repositioning sling can be pulled through. Ensure all loops are accessible and not under the patient's body.
- If rolling is not possible, use two flat slide sheets and panel these under the patient as per local manual handling policy. Once they are in place, the repositioning sling can be inserted between the slide sheets from head to foot.
- Once the repositioning sling is in place, the slide sheets can be removed from feet to head ensuring that all loops are accessible and not under the patient's body.
- Ensure that there are no creases under the patient's body. If there are, two caregivers should work at opposite corners of the sling and gently pull to remove creases.
- Ensure that loops do not come in contact with the floor. These can be tucked underneath the mattress if required.

REMOVING THE REPOSITIONING SLING

- Ensure the brakes of the bed/trolley are applied and the bed is of a suitable working height for the caregiver.
- The sling can be removed by turning the patient onto their side as per local manual handling policy, rolling the sling lengthways underneath itself and tucking under the patient's body.
- The patient can then be turned to the other side, the sling can be removed and patient returned to a supine position.
- If turning the patient is not possible, sling can also be removed by inserting two flat slide sheets.
- Slide sheets should be panelled above and below the repositioning sling and once inserted, the repositioning sling can be removed from feet to head.
- Once repositioning sling has been removed, slide sheets can then be removed. Follow your local manual handling policy for further guidance on paneling.

TURNING USING THE REPOSITIONING SLING

- Ensure the brakes of the bed/trolley are applied.
- Ensure the appropriate hoist and spreader bar have been selected.
- Position the spreader bar in line with the patient body and ensure the spreader bar is of suitable height to allow for loop attachment.
- If turning the patient to the left, select the right hand side loops. If turning the patient to the right, select the left hand side loops.
- Select appropriate inside or outside loop and attach to the spreader bar starting from the head loop. Colour loop selection can be recorded on the sling label or patients care plan.
- Ensure that the loops are correctly attached to the spreader bar before raising the hoist.
- Raise the hoist to turn the patient on to their side.
- Once in the correct position, a wedge or pillow can be placed at the patient's back for support.
- Once the patient is comfortable, lower the spreader bar and detach the loops.

MOVING UP THE BED

- Ensure the brakes of the bed/trolley are applied.
- Ensure the appropriate hoist and spreader bar have been selected.
- Attach either the inside or outside loops of the repositioning sling to the spreader bar starting from the head section.
- If possible, panel a slide sheet underneath the repositioning sling to reduce friction and shear.
- Ensure that the loops are correctly attached to the spreader bar before raising the hoist.
- Raise the hoist and move the patient up the bed to the desired position.
- Once the position has been achieved, lower the spreader bar, detach loops and remove slide sheet (if used).

LATERAL TRANSFER USING REPOSITIONING SLING

- Ensure the brakes of the bed/trolley are applied.
- Ensure the appropriate hoist and spreader bar have been selected.
- Attach either the inside or outside loops of the repositioning sling to the spreader bar starting from the head section.
- If not all loops can be attached, ensure that there is an equal number on each side to provide the best support to the patient.
- Ensure that loops are correctly attached to the spreader bar before raising hoist.
- Raise the hoist enough to lift the patient body clear off the surface ensuring that all wires, tubes, catheter bags etc are free to be moved.
- Patients can be transferred over onto the new surface ensuring that the brakes have been applied on the surface and the transfer gap between surfaces is small.
- Once the position has been achieved, lower the spreader bar, detach loops and remove slide sheet (if used).
- Lower the hoist to position the patient onto the new surface.
- Ensure that the patient is in the correct position and is comfortable.
- Lower the spreader bar, detach loops and remove hoist.

PRONE POSITION USING REPOSITIONING SLING

- Ensure the brakes of the bed/trolley are applied.
- Ensure the appropriate hoist and spreader bar have been selected.
- Insert a slide sheet underneath the repositioning sling using the panelling method.
- Insert a second repositioning sling under the slide sheet to ensure that the patient is prone onto a new sling.
- Position the spreader bar in line with the patient's body and ensure the spreader bar is of a suitable height to allow for loop attachment.
- If turning the patient to the left, select the right-hand side loops. If turning the patient to the right, select the left-hand side loops. Select the appropriate inside or outside loop and attach to the spreader bar starting from the head loop. Colour loop selection can be recorded on the sling label or the patient's care plan.
- Ensure that loops are correctly attached to the spreader bar before raising the hoist.
- Raise the hoist to turn the patient onto their side.
- Ensure the patient is comfortable before continuing
- Raise the hoist further whilst the caregivers on the receiving side support the patient into the prone position.
- Once the patient is in the desired position, the top repositioning sling and slide sheet can be removed.

RETURNING TO SUPINE POSITION FROM PRONE

- Ensure the brakes of the bed/trolley are applied.
- Ensure the appropriate hoist and spreader bar have been selected.
- Insert a slide sheet underneath the repositioning sling using the panelling method.
- Insert a second repositioning sling under the slide sheet to ensure that the patient is turned onto a new sling.
- Position the spreader bar in line with the patient's body and ensure the spreader bar is of a suitable height to allow for loop attachment.
- If turning the patient to the left, select the right-hand side loops. If turning the patient to the right, select the left-hand side loops.
- Select appropriate inside or outside loop and attach to the spreader bar starting from the head loop. Colour loop selection can be recorded on the sling label or the patient's care plan.
- Ensure that loops are correctly attached to the spreader bar before raising the hoist.
- Raise hoist to turn the patient onto their side.
- Ensure the patient is comfortable before continuing
- Raise the hoist further whilst the caregivers on the receiving side support the patient into the supine position.
- Once the patient is in the desired position, the top repositioning sling and slide sheet can be removed.

STORAGE AND HANDLING

Keep dry and out of direct sunlight.

DISPOSAL PROCEDURES

Follow local guidelines.

LIFETIME OF DEVICE

LOLER regulations require you to carry out a visual inspection of the sling prior to each use. The expected life of the repositioning sling is dependent on the actual conditions of use. If the sling shows any sign of damage or wear, dispose of it and select a new sling. If the sling is in contact with broken skin, do no use for longer than 30 days.

CLEANING

GBUK Banana Repositioning Slings are designed to be reused with the same patient i.e. the sling is single patient use only. If lightly soiled, manually wipe down the sling as soon as possible using alcohol based cleaning solution (up to 70%) or normal disinfectant solutions until the sling is visibly clean. If the sling does not appear visibly clean and remains contaminated upon visual inspection, discard the sling in line with local hospital protocol.

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GBUK Group award-winning specialist companies are experts in the provision of enteral feeding, critical care and patient moving and handling products and technologies. GBUK Group companies are suppliers to every UK NHS Trust and to a rapidly expanding export customer base.

GBUK Banana is a trading name of **GBUK Group Ltd.**

Symbol glossary available at www.gbukbanana.com.

GBUK Group Ltd.
Woodland House,
Blackwood Hall Business Park,
North Duffield, Selby,
North Yorkshire, YO8 5DD, UK
EU: +44 (0) 1757 282 945
US: 18665 289 422
info@gbukentral.com

QualRep Services BV
Utrechtseweg 310 B42
6812 AR Arnhem,
The Netherlands

B-REPSLING-03 2210 02