

# Comfort Value™ Stand Assist Slings

## INSTRUCTIONS FOR USE

### PRODUCT CODES & DESCRIPTIONS

All COMFORT VALUE™ Stand Assist Slings, product codes prefixed with:  
• CVSA - COMFORT VALUE™ Stand Assist

### INSTRUCTIONS FOR USE

Please read these instructions before use.

### DEVICE DESCRIPTION

The COMFORT VALUE™ Stand Assist slings are non-sterile medical devices, compliant with BS EN ISO 10535. Each sling has a printed label containing all pertinent information.

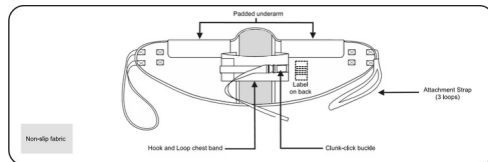


Figure 1

For further information, support and product training requests visit:  
[www.gbukgroup.com](http://www.gbukgroup.com)

### PATIENT INFORMATION

The following information is provided for use by trained healthcare professionals, and carers, and those with the appropriate training or knowledge, however as the learned intermediary between the company and the patient, the trained healthcare professionals and carers, and those with the appropriate training or knowledge must convey the aspects they consider relevant to the individual patient. The patient must be informed of the potential adverse effects & residual risks contained within this document.

### PATIENT TARGET GROUP(S)

Intended for patients who require lifting and/or repositioning following a suitable risk assessment.

### INTENDED USE

Intended to be used to transfer patients following a risk assessment. The sling should be used in conjunction with a mechanical lifting device (either mobile or ceiling hoist) to ensure safe lifting and repositioning of a patient.

### INTENDED USER(S)

The device is intended to be used by trained healthcare professionals and carers, and those with the appropriate training or knowledge. Refer to the Royal College of Nursing Guidelines.

### INDICATION

The sling is indicated for use on patients with limited mobility requiring lifting and/or repositioning.

### SAFE WORKING LOAD (SWL)

The SWL is provided on the label(s) directly attached to the device.

### COMPATIBILITY

Please visit [www.gbukgroup.com](http://www.gbukgroup.com) for a list of compatible hoists.

### CONTRAINDICATIONS

- Do not use on patients whose weight exceeds the safe working load limit of the sling or hoist, whichever is lower.

### WARNINGS & PRECAUTIONS

- A risk assessment must be carried out by a suitably qualified individual before any sling is selected and used, to ensure safety for the patient and handler. See 'RISK ASSESSMENT' section.
- Always carry out a visual check of the sling prior to every use for any damage, rips, tears, wearing, loose stitching, material degradation, or any other general signs of wear and tear. The sling should be disposed of in accordance with the disposal instructions, if any of these are present, and/or in the instance of failed LOLER inspection, and a new sling selected.
- Do not use with incompatible hoists.
- Follow local manual handling policy from your organization, National Back Exchange or HSE guidelines.
- All users must be trained or have the appropriate training or knowledge before using the sling.
- Inspect the device, packaging and labelling prior to use and do not use if damaged.
- Do not exceed SWL of the sling.
- The sling is subject to Lifting Operations and Lifting Equipment Regulations (LOLER) Examinations. A LOLER examination must be performed every 6 months, by trained and qualified persons.
- Wash the sling in accordance with the instructions provided in the 'CLEANING' section.
- Do not bleach, iron or dry clean the sling. Do not use stain removers.
- Do not exceed maximum washing and drying temperatures.
- Do not wash below minimum washing temperature.
- Do not place the sling on hot pipes or heaters. Keep away from flames.
- Where the sling has hooks and loops, ensure any hooks and loops are closed prior to laundering.
- Repeated washing may lead to sling deterioration and/or discolouration to some materials or fading of labels.

### EXPECTED LIFETIME

The expected performance and safety lifetime of the device is 5 years, however slings may last longer. All slings must be LOLER inspected and the inspection must be within date prior to use. Carry out a risk assessment prior to each use of the sling - Refer to 'WARNINGS & PRECAUTIONS' & 'RISK ASSESSMENT'.

### PERFORMANCE CHARACTERISTICS

- SWL as stated on label(s)
- BS EN ISO 21856 and BS EN ISO 10535 compliant

### DIRECTIONS FOR USE

Refer to 'Warnings & Precautions' and 'Risk Assessment' Sections prior to use of the sling.

- Ensure the patient is in a seated position, and a risk assessment and LOLER examination has been completed.
- With the label facing outwards, slide the sling down to midback, keeping it centrally aligned and with the users arms positioned to the underarm sections.  
NOTE: It may be necessary to support the user to prevent them tipping forward.
- With users arms on the outside, wrap the hook and loop chest band around the torso and fasten. Fasten the clunk-click buckle (there will be an audible click) and perform a 'tug test'. Tighten straps.
- Bring in the standing aid and adjust position as required and where model allows. Attach the appropriate strap loops onto the standing aid handle bar.
- Perform a 'tug test' and check for snagging. Slowly raise user into a standing position. To lower a user back into chair, repeat above steps 2-4 in reverse order.

\*A manual 'Tug Test' is recommended to ensure buckles and/or straps are properly secured prior to hoisting. For buckles, tug straps in opposite directions. For spreader bar attachment, tug straps in a downward motion.

### POTENTIAL ADVERSE EFFECTS & RESIDUAL RISKS

As with any procedure involving manual handling, use of manual handling aids, and hoisting equipment the typical risks associated apply. Failure of the device during use may expose the patient and user to risks such as falling, and potential subsequent associated injury, and/or user exposure to sudden mechanical stress and potential associated pain and/or musculoskeletal injuries.

### RISK ASSESSMENT

A risk assessment and visual inspection must be performed by competent person(s) before using the sling. This should consider the following factors:

- The physical and medical status of the patient, including skin integrity.
- The status of the sling, including the presence of any damage, rips, tears, soil, loose stitching, or any other general signs of wear and tear.
- The compatibility of the hoist and sling, including if the correct size and type of the sling is being used.
- The capability of the users, using the equipment; and,
- The situation in which it is being used

### WARRANTY

COMFORT VALUE™ slings are covered by a 3-year warranty. Any manufacturing defect will be rectified free of charge. Product damage through wear and tear, neglect or accidental damage is not covered.

### CLEANING & DISINFECTION

Routinely clean the sling, in accordance with these instructions and your facility's infection control protocols. Always follow your organisation's infection control guidelines to determine the appropriate cleaning schedule. Provide a copy of these instructions to the person(s) laundering the sling.

- Where present, close hook & loop fastenings before washing to prevent damage to the fabric and webbing.
- Wash gently at a minimum of 40°C and a maximum of 95°C. Use a mild soap solution or a non-biological detergent without bleaching additives.
- Do not add bleach or stain removers.
- Air dry, cabinet dry or tumble dry on a cool cycle (below 60°C). Do not tumble dry on a hot setting.

### STORAGE CONDITIONS

Standard storage conditions. Keep off the floor. Keep out of direct sunlight. Keep dry.

### DISPOSAL

No specific disposal requirements, handle contaminated items as clinical waste. Dispose of device in line with local and/or hospital policy.

### COMPLAINTS

Any user who has a complaint or experiences dissatisfaction in the product quality, identity, durability, reliability, safety, usability, effectiveness, and/or performance, should notify the manufacturer and distributor immediately.

In the event of a serious incident involving the device, the user should report the event to the manufacturer (GBUK Group Ltd) and the competent authority of the member state in which the user and/or patient is established, immediately, by telephone, email ([complaints@gbukgroup.com](mailto:complaints@gbukgroup.com)) or written correspondence.

When filing a complaint, provide the component(s) REF and SN number(s), your name and contact details, and the nature of the complaint.

In the event of a suspected device fault the user should report the event to GBUK Group Ltd with the above information and return the device to GBUK Group Ltd for examination.

COMFORT VALUE™ is a trade mark of Care & Independence Ltd., registered in the UK and other countries and regions.



#### GBUK Group Ltd.

Woodland House,  
Blackwood Hall Business Park,  
North Duffield, Selby,  
North Yorkshire,  
YO8 5DD, UK.

**T: EU: +44 (0) 1757 288587**

**E: [info@gbukgroup.com](mailto:info@gbukgroup.com)**

**W: [www.gbukgroup.com](http://www.gbukgroup.com)**



#### GBUK GmbH

Wirtschaftsprüfungsgesellschaft  
Im Klapperhof 7-23  
50670 Köln, Germany



#### QualRep Services BV

Utrechtseweg 310 B42,  
6812 AR, Arnhem,  
The Netherlands



#### Swiss AR Services AG

Industriestrasse 47,  
6300 Zug, Switzerland

